RIKTA®

KOH-3

KON-3

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CONTENTS

RU	3
EN	
INTENDED USE.....	12
GENERAL INFORMATION	12
SAFETY PRECAUTIONS	13
DELIVERY SET	14
TECHNICAL SPECIFICATIONS.....	14
OPERATING INSTRUCTION.....	16
NOZZLES DISINFECTING	17
PACKAGE, STORAGE AND TRANSPORTATION	17
MANUFACTURER'S WARRANTY.....	18
FR	19
DE	27
ACCEPTANCE CERTIFICATE	35
GUARANTEE CARD	35

KON-3 set of light-guide nozzles is certified in the Russian Federation jointly with RIKTA® device. Certificate of conformity № РОСС RU.ИМ02.Н17822.

INTENDED USE

The light-guide nozzles are designed to extend the functionality of the magnetic-infrared laser therapeutic devices of RIKTA® series through the ensuring of treatment zone by radiated light output.

EN

GENERAL INFORMATION

- KON-3 set of light-guide nozzles (hereinafter referred to as “the light-guide nozzles”) is a set of specialized light guide, which are used jointly with the magnetic-infrared laser therapeutic devices of RIKTA® series.
- The manufacturer reserves the right to make adjustments and amendments to this Manual and technical documentation at any time without notice.
- The manufacturer reserves the right to make changes in the configuration and design of the product, which will not affect its operation and basic technical specifications.

SAFETY PRECAUTIONS

- When using the light-guide nozzle to carry out a procedure of magnetic-infrared laser therapy a medical personnel and a patient is recommended to put on the special goggles.
- It is necessary to put a condom on the nozzle (to protect against infection) during the intracavitary exposure.
- **STRICTLY FORBIDDEN** to set/change the nozzles during the procedures!
- **STRICTLY FORBIDDEN** to direct emitting ends of the light-guide nozzles toward eyes, especially when irradiating the face.
- Chips, cracks on the light-guide nozzles which deteriorate its optical characteristics are not permissible.

EN

The nozzles disposal shall meet the conventional standards and local regulations for the environment protection.

DELIVERY SET


Description	Quantity
KON-3 set of light-guide nozzles:	
Light-guide nozzle "L"	1 pc
Light-guide nozzle "K"	1 pc
Light-guide nozzle "S"	1 pc
Operating manual	1 pc
Methodical manual	1 pc
Package	1 pc

EN

TECHNICAL SPECIFICATIONS

- The light-guide nozzles are made of non-toxic materials. The nozzles' surfaces are polished to reduce lateral radiation.
- The overall dimensions of nozzles (mm) and direction of light emission are shown on *Fig. 1*.
- Transmission factors of the light-guide nozzles are specified in Table:

Type of nozzle	Transmission factor, <i>no less</i>
Light-guide nozzle "L"	0,45
Light-guide nozzle "K"	0,85
Light-guide nozzle "S"	0,6

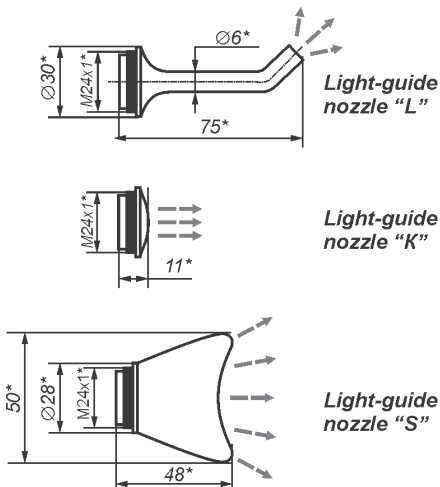EN

Fig. 1. KON-3 set of light-guide nozzles:

- Light-guide nozzle "L"** — used in gynecology and proctology;
- Light-guide nozzle "K"** — used in otorhinolaryngology and dentistry;
- Light-guide nozzle "S"** — used for laser puncture.

OPERATING INSTRUCTION

- Take the light-guide nozzle required for specific treatment procedure and screw it against stop into a threaded hole in the emitter aperture (the device must be switched off).

ATTENTION!

STRICTLY FORBIDDEN

mounting or replacement the nozzles during the procedures!

EN

ATTENTION!

Don't tighten the nozzle excessively to avoid the thread damage!

- The treatment procedures shall be carried out in accordance with the "Methodical recommendations on application".

ATTENTION!

It is necessary to put a condom on the nozzle during **the intracavitary exposure**

NOZZLES DISINFECTING

ATTENTION!

STRICTLY FORBIDDEN

to use alcohol-containing means
for disinfection of nozzles!

- **When using the light-guide nozzles at home:** it is recommended to wipe the nozzle using a **non-alcoholic antiseptic** soft cloth or gauze /cotton wetted in a special non-alcoholic cleanser for disinfecting of medical appliances.
- **When using the light-guide nozzles in health care facilities:** it is recommended preliminary cleaning of the nozzle at room temperature during 3 hours in 0.5% peroxide solution by adding 0.5% cleanser.

EN

PACKAGE, STORAGE AND TRANSPORTATION

- The light-guide nozzles shall be packed in the Manufacturer's package.
- The light-guide nozzles should be stored in the Manufacturer's package in the heated room at temperature in the range of +10...+35 °C and relative humidity of no more than 80% (at +25 °C).
- The light-guide nozzles should be transported by all types of the closed transport in Manufacturer's package. It is necessary to ensure a stable position of the package during transportation.

MANUFACTURER'S WARRANTY

- The Manufacturer guarantees that the KON-3 set of light-guide nozzles parameters correspond to the requirements of the Technical Specifications TU 9444-007-18537172-2003, within 24 months from the date of sale (if the consumer complies with the conditions of use, transportation and storage).
- In the absence of the sales date and the seller's stamp in the Sale certificate, the warranty period shall be calculated from the date of Manufacturing of light-guide nozzles.
- All issues related to repairs and maintenance should be directed to the Manufacturer or to the technical centers listed on the website www.rikta.ru or in a separate warranty card (if available).

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